




Air Filter

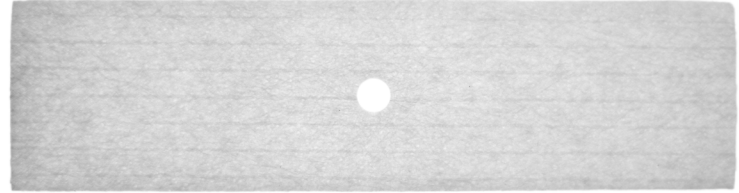

For use with the Giraffe™ Incubator Carestation and the Giraffe™ Omnibed Carestation.

Maintain incubator airflow while providing high-efficiency filtration

Features

- Simple installation
- Bi-directional, fits in the cover in any direction
- Blank label to mark installation date and track for replacement



Air Filter, bi-directional

**Our equipment.
Our accessories.
Simply compatible.**

You can trust GE Healthcare accessories and consumables to keep your clinical devices working efficiently.

- Simplified ordering process
- Fast delivery
- One contact point for expert advice
- Quality support at every step

PRODUCT SPECIFICATIONS

	Part Information
Micro Filter	0.3 μ -99.8% efficiency @ flow rate of 1 cm/sec
System Performance	CO ₂ level \leq 0.3%
Weight	3.54 to 5.31 oz. (120 to 180 GSM)
Dimensions	5.58" x 1.43" x 0.08" (L x W x T) approximate

Materials

Electrostatic synthetic, non-woven

100% synthetic fibers

PP spunbond, non-woven

Product use

Air filter should be changed every 3 months. In addition, it should be changed after the bed was used by an infected patient or when visibly dirty.

Label the filter site with the date of last change and inspect filter at every cleaning.

Operating environment

Temperature:	20° to 30°C (68° to 86°F)
Humidity:	\leq 85% RH (non-condensing)
Pressure:	50-106kPa

Storage

Temperature:	-25° to 60°C (-13° to 140°F)
Humidity:	\leq 85% RH (non-condensing)
Pressure:	50-106kPa

Product may not be available in all countries and regions. Contact a GE Healthcare Representative for more information. Please visit www.gehealthcare.com.

Imagination at work

© 2020 General Electric Company – All rights reserved.

GE Healthcare reserves the right to make changes in specifications and features shown herein, or discontinue the product described at any time without notice or obligation. Contact your GE Healthcare representative for the most current information. GE and the GE Monogram, are trademarks of General Electric Company. GE Healthcare, a division of General Electric Company. BiliSoft, BiliBlanket and Giraffe are trademarks of General Electric Company. GE Medical Systems, Inc., doing business as GE Healthcare. All other product names and logos are trademarks or registered trademarks of their respective companies.

DOC2390726 Rev 1 5/20

Disposal considerations

Dispose as clinical waste, in accordance with hospital policy, local guidelines and regulations.

Packaging materials

Primary: Polybag – Polyethylene (PE)

Compliance

Not a medical device.


Regulatory standards

Air filter compliant with REACH regulation (EC)/1907/2006.

Air filter compliant with (RoHS) directive, 2011/65/EU, (Ref: GEHC ROHS compliant per 5240305GSP).

PART LIST

Part Number	Description	Quantity
2107356-001	Air filter, bi-directional	10

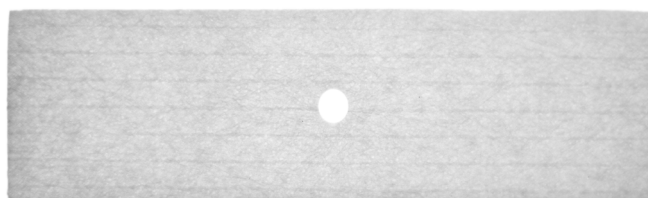

© 2020 General Electric Company – All rights reserved.

GE Healthcare reserves the right to make changes in specifications and features shown herein, or discontinue the product described at any time without notice or obligation. Contact your GE Healthcare representative for the most current information. GE and the GE Monogram, are trademarks of General Electric Company. GE Healthcare, a division of General Electric Company. BiliSoft, BiliBlanket and Giraffe are trademarks of General Electric Company. GE Medical Systems, Inc., doing business as GE Healthcare. All other product names and logos are trademarks or registered trademarks of their respective companies.

DOC2390726 Rev 1 5/20