




Infant Resuscitator T-Piece Circuit

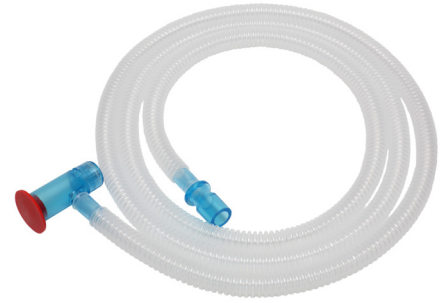

Disposable resuscitator T-piece circuits with adjustable PEEP

Disposable Resuscitator T-Piece Circuits, designed for neonates, provide accurate, reliable ventilation and pulmonary resuscitation. They are suitable for labor wards, delivery suites, NICU, SCBU, and transport.

Features

- The adjustable positive end expiratory pressure (PEEP) control valve easily rotates to maintain the exact pressure setting. Pressure remains the same even during extended use
- A comfortable, ergonomic design can help reduce hand fatigue when used for extended periods
- The blue-colored T-piece makes it easily visible in delivery suites, in the NICU and during transport of patients
- A thumb/finger port enables controlled inspiration and expiration. Flow can be controlled to deliver 21% to 100% oxygen or free flow by simply placing the thumb over the port
- Made from lightweight, flexible and crush resistant tubing, it can maintain the flow of air even when twisted

Product compatibility

With a universal connector, it fits with all commonly used manual resuscitator units and does not require an adapter.

Our equipment. Our accessories. Simply compatible.

You can trust GE Healthcare accessories and consumables to keep your clinical devices working efficiently.

- Simplified ordering process
- Fast delivery
- One contact point for expert advice
- Quality support at every step

PRODUCT SPECIFICATIONS

Materials

Component	Material
PEEP valve cap	Acrylonitrile Butadiene Styrene (ABS)
O ring	Silicone
Elbow	Acrylonitrile Butadiene Styrene (ABS)
15mm Male manifold cap	Polypropylene (PP)
10mm Corrugated tubing	Low Density Polyethylene (LDPE)
10M-15M/15F Connector	Acrylonitrile Butadiene Styrene (ABS) Does not contain natural rubber latex, or DEHP

Operating conditions

Check the host device manual for recommendations.

PART LIST

Part Number	Description	Quantity
2109059-001	Disposable T-Piece Patient Circuit Kit	20
2109060-001	Disposable T-Piece Patient Circuit Kit with Mask Size 0	10
2109061-001	Disposable T-Piece Patient Circuit Kit with Mask Size 1	10

Storage conditions

Store in a cool, dry place out of direct sunlight.

Shelf life

The shelf life is five years from the date of manufacture. This is based on the stability of the devices' components and raw materials sourced. Expiration date is clearly marked on individual product pouch.

Compliance

Connectors compliant with ISO 5356

Compatible with ISO 5356 compliant resuscitation units, neonatal face masks, ET tubes.

Imagination at work

Product may not be available in all countries and regions. Contact a GE Healthcare Representative for more information. Please visit www.gehealthcare.com.

Distributed by:

Datex-Ohmeda, Inc.



Manufactured by:

Flexicare Medical Ltd.
Cynon Valley Business Park
Mountain Ash CF45 4ER
UK